

**Human Peroxisome Proliferator-Activated
Receptor Delta
(NR1C2, PPARD, PPAR δ , PPAR β)**

Reporter Assay System

3x 32 Assays in 96-well Format
Product # IB00121-32

▪


Technical Manual
(version 7.1c)

www.indigobiosciences.com

3006 Research Drive, Suite A1, State College, PA 16801, USA

Customer Service:
814-234-1919; FAX 814-272-0152
customerserv@indigobiosciences.com

Technical Service:
814-234-1919
techserv@indigobiosciences.com

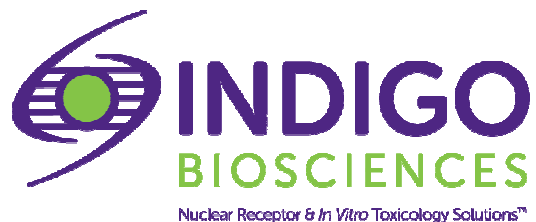

Human PPAR δ Reporter Assay System 3x 32 Assays in 96-well Format

I. Description	
▪ The Assay System.....	3
▪ The Assay Chemistry.....	3
▪ Preparation of Test Compounds.....	4
▪ Assay Scheme.....	4
▪ Assay Performance.....	5
II. Product Components & Storage Conditions	6
III. Materials to be Supplied by the User	6
IV. Assay Protocol	
▪ A word about <i>Antagonist</i> -mode assay setup.....	7
▪ <i>DAY 1 Assay Protocol</i>	7
▪ <i>DAY 2 Assay Protocol</i>	9
V. Related Products	10
VI. Limited Use Disclosures	11
APPENDIX 1: Example Scheme for Serial Dilution	12

I. Description

▪ The Assay System ▪

This nuclear receptor assay utilizes proprietary non-human cells engineered to provide constitutive, high-level expression of the **Human Peroxisome Proliferator-Activated Receptor Delta** (NR1C2), a ligand-dependent transcription factor commonly referred to as PPAR δ or **PPAR δ** .

INDIGO's Reporter Cells include the luciferase gene functionally linked to a PPAR δ -responsive promoter. Thus, quantifying changes in luciferase expression in the treated cells provides a sensitive surrogate measure of the changes in PPAR δ activity. The principal application of this assay is in the screening of test samples to quantify any functional activity, either agonist or antagonist, that they may exert against human PPAR δ .

PPAR δ Reporter Cells are prepared using INDIGO's proprietary **CryoMite™** process. This cryo-preservation method yields exceptional cell viability post-thaw, and provides the convenience of immediately dispensing healthy, division-competent cells into assay plates. There is no need for cumbersome intermediate treatment steps such as spin-and-rinse of cells, viability determinations, cell titer adjustments, or the pre-incubation of cells prior to assay setup.

INDIGO Bioscience's Nuclear receptor assay kits provide all-inclusive cell-based assay systems. In addition to PPAR δ Reporter Cells, this kit provides two optimized media for use during cell culture and in diluting the user's test samples, a reference agonist, Luciferase Detection Reagent, and a cell culture-ready assay plate.

▪ The Assay Chemistry ▪

INDIGO's assays capitalize on the extremely low background, high-sensitivity, and broad linear dynamic range of bio-luminescence reporter gene technology.

Reporter Cells incorporate the cDNA encoding beetle luciferase, a 62 kD protein originating from the North American firefly (*Photinus pyralis*). Luciferase catalyzes the mono-oxidation of D-luciferin in a Mg⁺²-dependent reaction that consumes O₂ and ATP as co-substrates, and yields as products oxyluciferin, AMP, PP_i, CO₂, and photon emission. Luminescence intensity of the reaction is quantified using a luminometer and is reported in terms of Relative Light Units (RLU's).

INDIGO's Assays feature a luciferase detection reagent specially formulated to provide stable light emission between 5 and 90+ minutes after initiating the luciferase reaction. Incorporating a 5-minute reaction-rest period ensures that light emission profiles attain maximal stability, thereby allowing assay plates to be processed in batch. By doing so, the signal output from all sample wells, from one plate to the next, may be directly compared within an experimental set.


▪ Preparation of Test Compounds ▪

Most commonly, test compounds are solvated at high concentration in DMSO, and these are stored as master stocks. Master stocks are then diluted to appropriate working concentrations immediately prior to setting up the assay. Users are advised to dilute test compounds to 2x-concentration stocks using **Compound Screening Medium (CSM)**, as described in *Step 2* of the **Assay Protocol**. This method avoids the adverse effects of introducing high concentrations of DMSO into the assay. The final concentration of total DMSO carried over into assay reactions should never exceed 0.4%.

NOTE: CSM is formulated to help stabilize hydrophobic test compounds in the aqueous environment of the assay mixture. Nonetheless, high concentrations of extremely hydrophobic test compounds diluted in CSM may lack long-term stability and/or solubility, especially if further stored at low temperatures. Hence, it is recommended that test compound dilutions are prepared in CSM immediately prior to assay setup and are considered to be 'single-use' reagents.

▪ Assay Scheme ▪

Figure 1. Assay workflow. *In brief*, Reporter Cells are dispensed into wells of the assay plate and then immediately dosed with the user's test compounds. Following 22 -24 hr incubation, treatment media are discarded and prepared Luciferase Detection Reagent (LDR) is added. Light emission from each assay well is quantified using a plate-reading luminometer.

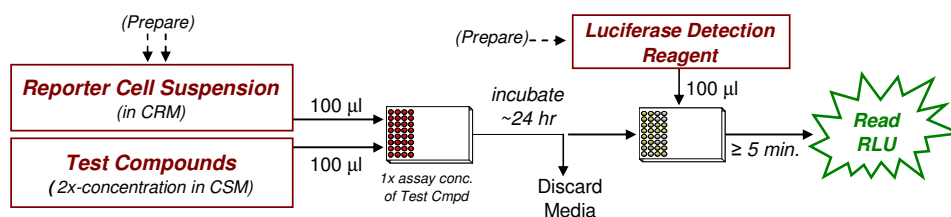

▪ Assay Performance ▪

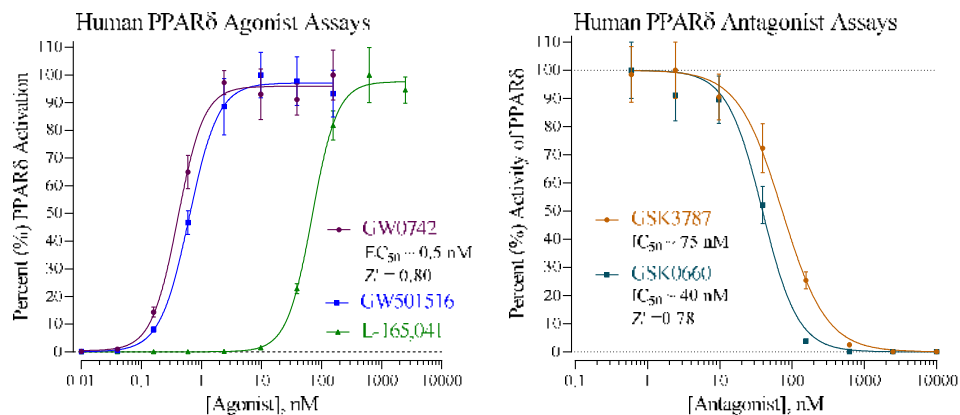


Figure 2. Agonist and Antagonist dose-response of the PPAR δ receptor.

PPAR δ reference agonists GW0742 (provided), GW501516 and L-165041 (Tocris), and antagonists GSK3787 and GSK0660 (Tocris) were analyzed. Average Relative Light Units (RLU) and their respective values of Standard Deviation (SD), Coefficient of Variation (CV), Fold-Activation and % Activity were calculated for each treatment concentration (n = 4). Z' values were calculated as per Zhang, *et al.* (1999)¹.

Data plots were generated using GraphPad Prism software. Compound concentrations were Log transformed and RLU values were normalized to percent (%) activity of the receptor. Non-linear regression analyses were performed and EC₅₀ & IC₅₀ values determined.

The large dynamic response and Z' values confirm the robust performance of the PPAR δ agonist and antagonist assays and demonstrate their suitability for use in HTS applications.

¹ Zhang JH, Chung TD, Oldenburg KR. (1999) A Simple Statistical Parameter for Use in Evaluation and Validation of High Throughput Screening Assays. *J Biomol Screen.*:4(2), 67-73.

$$Z' = 1 - [3 * (SD^{\text{Reference}} + SD^{\text{Vehicle Bkg}}) / (RLU^{\text{Reference}} - RLU^{\text{Vehicle Bkg}})]$$

II. Product Components & Storage Conditions

This Human PPAR δ Assay kit contains materials to perform three distinct groups of assays in a 96-well plate format. Reagents are configured so that each group will comprise 32 assays. If desired, however, reagents may be combined to perform either 64 or 96 assays.

The individual aliquots of Reporter Cells are provided as single-use reagents. Once thawed, cells can NOT be refrozen or maintained in extended culture with any hope of retaining downstream assay performance. Therefore, extra volumes of these reagents should be discarded after assay setup.

Assay kits are shipped on dry ice. Upon receipt, individual kit components may be stored at the temperatures indicated on their respective labels. Alternatively, the entire kit may be further stored at -80°C.

To ensure maximal viability, "Reporter Cells" must be maintained at -80°C until immediately prior to use.

The date of product expiration is printed on the Product Qualification Insert (PQI) enclosed with each kit.

<u>Kit Components</u>	<u>Amount</u>	<u>Storage Temp.</u>
▪ PPAR δ Reporter Cells	3 x 0.60 mL	-80°C
▪ Cell Recovery Medium (CRM)	1 x 10.5 mL	-20°C
▪ Compound Screening Medium (CSM)	1 x 35 mL	-20°C
▪ GW0742, 1.0 mM (in DMSO) (reference agonist for PPAR δ)	1 x 30 μ L	-20°C
▪ Detection Substrate	3 x 2.0 mL	-80°C
▪ Detection Buffer	3 x 2.0 mL	-20°C
▪ Plate frame	1	ambient
▪ Snap-in, 8-well strips (white, sterile, cell-culture ready)	12	ambient

III. Materials to be Supplied by the User

The following materials must be provided by the user, and should be made ready prior to initiating the assay procedure:

DAY 1

- dry ice bucket
- cell culture-rated laminar flow hood.
- 37°C, humidified 5% CO₂ incubator for mammalian cell culture.
- 37°C water bath.
- 70% alcohol wipes
- 8- or 12-channel electronic, repeat-dispensing pipettes & sterile tips
- disposable media basins, sterile.
- sterile multi-channel media basins (such as the Heathrow Scientific "Dual-Function Solution Basin"), or deep-well plates, or appropriate similar vessel for generating dilution series of reference compound(s) and test compound(s).
- *Optional*: antagonist reference compound.
- *Optional*: clear 96-well assay plate, sterile, cell culture treated, for viewing cells on Day 2.

DAY 2 plate-reading luminometer.

IV. Assay Protocol

Review the entire Assay Protocol before starting. Completing the assay requires an overnight incubation. *Steps 1-8* are performed on **Day 1**, requiring less than 2 hours to complete. *Steps 9-14* are performed on **Day 2** and require less than 1 hour to complete.

▪ A word about Antagonist-mode assay setup ▪

Receptor inhibition assays expose the Reporter Cells to a constant, sub-maximal concentration (typically between EC₅₀ – EC₈₅) of a known agonist AND the test compound(s) to be evaluated for antagonist activity. This PPAR δ Assay kit includes a 1.0 mM stock solution of **GW0742**, an agonist of PPAR δ that may be used to setup antagonist-mode assays. 4 nM GW0742 typically approximates EC₈₀ in this assay. Hence, it presents a reasonable assay concentration of agonist to be used when screening test compounds for inhibitory activity.

We find that adding the reference agonist to the bulk suspension of cells (*i.e.*, prior to dispensing into assay wells) is the most efficient and precise method of setting up antagonist assays, and it is the method presented in *Step 5b* of the following protocol. Note that, in *Step 6*, 100 μ l of treatment media is combined with 100 μ l of pre-dispensed [Cells + agonist]. Consequently, one must prepare the bulk suspension of Reporter Cells to contain a 2x-concentration of the reference agonist. **APPENDIX 1** provides a dilution scheme that may be used as a guide when preparing cell suspension supplemented with a desired 2x-concentration of agonist.

DAY 1 Assay Protocol: All steps must be performed using aseptic technique.

1.) Remove Cell Recovery Medium (CRM) and Compound Screening Medium (CSM) from freezer storage and thaw in a 37°C water bath.

2.) Prepare dilutions of treatment compounds (first see *Note 5.3*): Prepare Test Compound treatment media for *Agonist-* or *Antagonist-mode* screens.

Total DMSO carried over into assay reactions should never exceed 0.4%.

Note that, in *Step 6*, 100 μ l of the prepared treatment media is added into assay wells that have been pre-dispensed with 100 μ l of Reporter Cells. Hence, to achieve the desired *final* assay concentrations one must prepare treatment media with a 2x-concentration of the test and reference material(s). Use **CSM** to prepare the appropriate dilution series. Manage dilution volumes carefully. This assay kit provides 35 ml of CSM.

Preparing the positive control: This PPAR δ Assay kit includes a 1.0 mM stock solution of the reference agonist **GW0742**. The following 8-point treatment series, with concentrations presented in 3-fold decrements, provides a suitable dose-response: 90, 30, 10, 3.33, 1.11, 0.370, 0.123, and 0.0411 nM, and including a 'no treatment' control. **APPENDIX 1** provides an example for generating such a dilution series.

3.) Rapid Thaw of the Reporter Cells: *First*, retrieve the tube of **CRM** from the 37°C water bath and sanitize the outside with a 70% ethanol swab.

Second, retrieve **Reporter Cells** from -80°C storage and immerse in dry ice to transport the tube to a laminar flow hood: 1 tube for 32 assay wells, 2 tubes for 64 assay wells, and 3 tubes for 96 assay wells. Transfer the tube(s) into a rack and, *without delay*, perform a rapid thaw of the frozen cells by transferring a 3.0 ml volume of 37°C CRM into each tube of frozen cells. Recap the tube of cells and immediately place it in a 37°C water bath for 5 - 10 minutes. If only one tube of reporter cells is thawed (32 assays), the resulting volume of cell suspension will be 3.6 ml.

(continued ...)

Third, work in the cell culture hood to *carefully* mount four sterile 8-well strips into the blank assay plate frame. Strip-wells are fragile. Note that they have keyed ends (square and round), hence, they will fit into the plate frame in only one orientation.

4.) Retrieve the tube of cell suspension from the water bath and sanitize the outside surface of the tube with a 70% alcohol swab.

5.) *a. Agonist-mode assays.* Gently invert the tube of cells several times to disperse cell aggregates and gain a homogenous cell suspension. Without delay, dispense 100 μ l of cell suspension into each well of the assay plate.

~ or ~

b. Antagonist-mode assays. Gently invert the tube of cells several times to disperse any cell aggregates, and to gain a homogenous cell suspension. Supplement the bulk suspension of cells with the desired 2x-concentration of reference agonist (refer to "A word about antagonist-mode assay setup", pg. 7). Dispense 100 μ l of cell suspension into each well of the assay plate.

NOTE 5.1: Take special care to prevent cells from settling during the dispensing period. Allowing cells to settle during the transfer process, and/or lack of precision in dispensing uniform volumes across the assay plate *will* cause well-to-well variation (= increased Standard Deviation) in the assay.

NOTE 5.2: Users sometimes wish to examine the cells using a microscope. If so, the extra volume of cell suspension provided with each kit may be dispensed (100 μ l/well) into a clear 96-well cell culture treated assay plate, followed by 100 μ l/well of CSM. Incubated overnight in identical manner to those cells contained in the white assay plate.

NOTE 5.3: For logistical reasons, some users find it more convenient to first plate the cells and then prepare their test compound dilutions. That strategy works equally well. Once plated, cells may be placed in an incubator for up to 3 hours before proceeding to *Step 6*.

6.) Dispense 100 μ l of 2x-concentration treatment media into appropriate assay wells.

7.) Transfer the assay plate into a 37°C, humidified 5% CO₂ incubator for 22 - 24 hours.

NOTE: Ensure a high-humidity (\geq 85%) environment within the cell culture incubator. This is critical to prevent the onset of deleterious "edge-effects" in the assay plate.

8.) For greater convenience on Day 2, retrieve the appropriate number of vials of **Detection Substrate** and **Detection Buffer** from freezer storage and place them in a dark refrigerator (4°C) to thaw overnight.

DAY 2 Assay Protocol: Subsequent manipulations do *not* require special regard for aseptic technique and may be performed on a bench top.

9.) 30 minutes before intending to quantify PPAR δ activity, remove **Detection Substrate** from the refrigerator and place them in a low-light area so that it may equilibrate to room temperature. Gently invert the tube several times to ensure a homogenous solution.

NOTE: Do NOT actively warm Detection Substrate above room temperature. If these solutions were not allowed to thaw overnight at 4°C, a room temperature water bath may be used to expedite thawing.

10.) Set the plate-reader to "luminescence" mode. Set the instrument to perform a single 5 second "plate shake" prior to reading the first assay well. Read time may be set to 0.5 second (500 mSec) per well, *or less*.

11.) *Immediately before proceeding to Step 12:* To read 32 assay wells, transfer the entire volume of 1 vial of Detection Buffer into 1 vial of Detection Substrate, thereby generating a 4 ml volume of **Luciferase Detection Reagent (LDR)**. Mix gently to avoid foaming.

12.) After 22-24 hours of incubation, remove media contents from each well.

NOTE: Because the assay plate is composed of a frame with snap-in strip-wells, the practice of physically ejecting media is NOT advised. Do *not* touch the well bottom, or run the tip of the aspiration device around the bottom circumference of the assay well. Such practices will result in destruction of the cells and greatly increased well-to-well variability. Complete removal of the media is efficiently performed by tilting the plate on edge and aspirating media using an 8-pin manifold (*e.g.*, Wheaton Science Microtest Syringe Manifold, # 851381) affixed to a vacuum-trap apparatus.

13.) Add 100 μ l of **LDR** to each well of the assay plate. Allow the assay plate to rest at room temperature for at least 5 minutes. Do not shake the assay plate during this period.

14.) Quantify luminescence.

V. Related Products

<i>Product No.</i>	<i>Product Descriptions</i>
Cynomolgus Monkey PPARδ Assay Products	
C00121-32	Monkey PPAR δ Reporter Assay 3x 32 assays in 8-well strips (96-well plate format)
C00121	Monkey PPAR δ Reporter Assay 1x 96-well format assay
Dog PPARδ Assay Products	
D00121-32	Dog PPAR δ Reporter Assay 3x 32 assays in 8-well strips (96-well plate format)
D00121	Dog PPAR δ Reporter Assay 1x 96-well format assay
Rat PPARδ Assay Products	
R00121-32	Rat PPAR δ Reporter Assay 3x 32 assays in 8-well strips (96-well plate format)
R00121	Rat PPAR δ Reporter Assay 1x 96-well format assay
Mouse PPARδ Assay Products	
M00121-32	Mouse PPAR δ Reporter Assay 3x 32 assays in 8-well strips (96-well plate format)
M00121	Mouse PPAR δ Reporter Assay 1x 96-well format assay
Human PPARδ Assay Products	
IB00121-32	Human PPAR δ Reporter Assay 3x 32 assays in 8-well strips (96-well plate format)
IB00121	Human PPAR δ Reporter Assay 1x 96-well format assay
IB00122	Human PPAR δ Reporter Assay 1x 384-well format assays
Panel of Human PPAR Assays	
IB00131-32P	PANEL of Human PPAR γ , PPAR α and PPAR δ Assays 32 assays each in 8-well strips (96-well plate format)

LIVE Cell Multiplex (LCM) Assay	
<i>Product No.</i>	<i>Product Descriptions</i>
LCM-01	Reagent volumes sufficient to perform 96 Live Cell Assays in 1x96-well, or 2x48-well, or 3x32-well assay plate formats
LCM-05	Reagent in 5x-bulk volume to perform 480 Live Cell Assays performed in 5x 96-well plates
LCM-10	Reagent in 10x-bulk volume to perform 960 Live Cell Assays performed in 10x 96-well plates

Please refer to INDIGO Biosciences website for updated product offerings.

www.indigobiosciences.com

VI. Limited Use Disclosures

Products commercialized by INDIGO Biosciences, Inc. are for RESEARCH PURPOSES ONLY – not for therapeutic, diagnostic, or contact use in humans or animals.

“CryoMite” is a Trademark TM of INDIGO Biosciences, Inc. (State College, PA, USA)

Product prices, availability, specifications and claims are subject to change without prior notice.

Copyright © INDIGO Biosciences, Inc. (State College, PA, USA). All Rights Reserved.

APPENDIX 1

Example scheme for the serial dilution of GW0742 reference agonist, and the setup of a PPAR δ dose-response assay.

